


MRI Transport Circuit

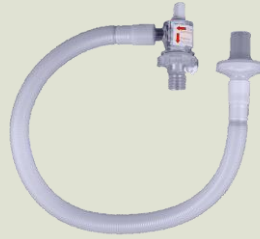

MRI Transport
Made Simple

Westmed
When Your Care Is Critical

MRI Transport Circuit



**Integral PEEP Valve,
4 – 20 cmH₂O**



**Also Available with CO₂ Easy™
CO2 Detector**



**Monitoring Port
with Tethered Cap**

A complete, lightweight, ready-to-use, non-metallic, non-electrically conductive MRI transport circuit, individually packaged. Compatible with Westmed BagEasy™ Manual Resuscitator and with some portable transport ventilators that have a passive exhalation valve. Features a minimal deadspace non-rebreathing valve with a swivel patient connector, adjustable PEEP valve and a proximal monitoring port that can be used for airway pressure. Available with Adult CO₂ Easy™ Carbon Dioxide Detector: 2-inch diameter indicator is easy to read and the anti-fog lens reduces condensation from exhaled gas. Purple-to-yellow color change allows exhaled CO₂ to be confirmed from a distance.

Benefits:

- Lightweight, decreases pull on ET tube
- MRI-Compatible
- Expandable tubing, ~28 – 96 inches
- Minimal Deadspace / Integral PEEP Valve
- Patient swivel connector, easy to attach
- Pre-assembled for simplicity in setup and use

Specifications:

- PEEP Valve Range (Calibrated Scale) = 4 to 20 cm H₂O (±4 cm H₂O); Minimum <3 cm H₂O; Maximum >20 cm H₂O.
- Deadspace of non-rebreathing valve = 10 mL.
- Tubing Length = ~28 inches, collapsed; ~8 feet (96 inches), extended.
- Minimum tubing compression/compliance (fully collapsed) = 0.47 mL/cm H₂O.
- Maximum tubing compression/compliance (fully extended) = 0.95 mL/cm H₂O.
- Bacterial/viral ventilator outlet filter included.



ORDERING INFORMATION

REF 9317 • MRI Ventilator Circuit, 96" MRI Disposable Circuit with Patient Valve and PEEP, 10/case

REF 9318 • MRI Ventilator Circuit, 96" MRI Disposable Circuit with Patient Valve, No PEEP, 10/case

Compatible with Westmed BagEasy™ Manual Resuscitator and with portable transport ventilators that have a passive exhalation valve.

Westmed

When Your Care Is Critical

Westmed, Inc.

5580 S. Nogales Highway

Tucson, AZ 85706

Toll Free: (800) 975-7987

Phone: (520) 294-7987

Fax: (520) 294-6061

Distributed by:



Please visit us at
www.westmedinc.com
for a complete listing of
Westmed products.

© 2017 Westmed, Inc. All rights reserved.
PN 76002, Rev. 04